

vetscan cPL RAPID TEST:

THE BEST CORRELATION TO THE **GOLD STANDARD** OF PANCREATITIS IN A PEER-REVIEWED STUDY, WITH THE CONVENIENCE OF IMMEDIATE AND **ACCURATE SEMI-QUANTITATIVE RESULTS**.

AN INDEPENDENT STUDY COMPARED THE TOP 4 CANINE PANCREATITIS ASSAYS¹


In-clinic
VETSCAN
cPL Rapid Test
&
SNAP[®] cPL[™]


Send-out
Spec cPL[®]
&
Precision PSL[™]

TO LOOK AT INTRACLASS CORRELATION COEFFICIENT (ICC)

Agreement Scoring

How strongly do these tests resemble each other?

On a scale from 0 to 1



Poor
<0.40



Fair
0.40 - 0.59



Good
0.60 - 0.74



Excellent
0.75 - 1.00



50 canines with gastrointestinal clinical signs

THEY TESTED



Collected serum samples



Tested serum on all 4 tests



SNAP cPL
Performed in-clinic



VETSCAN cPL Rapid Test
Performed in-clinic



Spec cPL
Sent out to Texas A&M University, Gastrointestinal Laboratory, College Station, TX.



Precision PSL
Sent out to Antech[®] Diagnostics, Irvine, CA.

AND FOUND



Poor
<0.40



Fair
0.40 - 0.59



Good
0.60 - 0.74



Excellent
0.75 - 1.00

All assays demonstrated **excellent** agreement with each other

ICC Score¹

0.86 SNAP cPL vs. Precision PSL

0.89 Spec cPL vs. Precision PSL

0.91 VETSCAN cPL vs. Precision PSL

0.92 SNAP cPL vs. Spec cPL

0.93 VETSCAN cPL vs. SNAP cPL

0.96 VETSCAN cPL vs. Spec cPL



BEST AGREEMENT WAS BETWEEN THE VETSCAN cPL RAPID TEST AND THE *GOLD STANDARD* OF PANCREATITIS ASSAYS, SPEC cPL

YOU CAN HAVE GREAT CONFIDENCE IN THE VETSCAN cPL RAPID TEST WHEN THERE IS A CLINICAL SUSPICION OF CANINE PANCREATITIS.*

Reference

¹ Cridge, H., MacLeod, AG., Pachtiger, GE., Mackin, A.J., Sullivan, AM., Thomason, JM., Archer, TM., Lunsford, KV., Rosenthal, K., Willis, RW. Evaluation of SNAP cPL, Spec cPL, VETSCAN cPL Rapid Test, and Precision PSL Assays for the Diagnosis of Clinical Pancreatitis in Dogs. J Vet Intern Med. 2018 Mar; 32(2):658-664. doi: 10.1111/jvim.15039. Epub 2018 Feb 9.

vetscan



COMPLETE POINT-OF-CARE WORKFLOW

The VETSCAN cPL Rapid Test is an easy-to-use, semi-quantitative test for the detection of cPL in serum or plasma at the point of care.


- + **Semi-quantitative results** - Convenient semi-quantitative immunoassay for the detection of canine pancreas-specific lipase in serum or plasma, that gives you easy-to-read numeric results.
- + **Results** - In a study, demonstrated excellent agreement with Spec cPL reference lab test.¹
- + **Automated interpretation** - Read results with the VETSCAN VUE.
- + **Room temperature storage** - Save refrigerator space and eliminate test warm-up time.¹
- + **Easy to use** - 3 simple steps help eliminate unnecessary re-tests and minimize the margin for human error.



"THE VETSCAN CPL IS NOT ONLY RAPID AND EASY TO USE, BUT GETTING AN ACTUAL NUMBER TAKES OUT THE GUESSING GAME AND MAKES PATIENT MONITORING AND FOLLOW-UP MUCH EASIER—GREATLY INCREASING MEDICAL QUALITY!"

- Dr. Philippe Roof, Pet Medical Center Boca

CONVENIENT. RAPID. EASY.

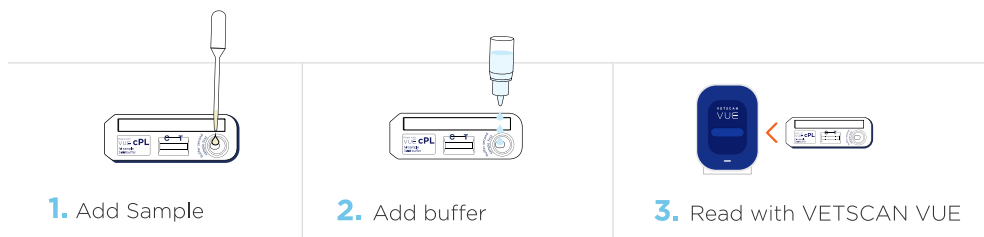

Contact Your Sales Representative Today!

888 963 8471

zoetisUS.com/vetscancpl

*Results of any of the 4 diagnostic assays alone, without supporting clinical findings are insufficient to establish a diagnosis of pancreatitis. IDEXX, SNAP, and SPEC cPL are registered trademarks of IDEXX Laboratories, Inc. Antech is a registered trademark of Antech Diagnostics. Precision PSL is a trademark of Antech Diagnostics. All trademarks are the property of Zoetis Services LLC or a related company or a licensor unless otherwise noted. ©2020 Zoetis Services LLC. All rights reserved. ABX-00122A R1

ZOETIS
PETCARE